


BALANCE ACCESSORIES

MOTORIZED CALIBRATION

Fully automatic calibration accomplished by a smooth noiseless motor driven standard weight inside the balance eliminating weight maintenance. Accurately machined stainless steel standard weight is used for the purpose. Any tolerance in the weight machining is corrected through software. Balance can be programmed to initiate the motorized calibration process every time when it is switched ON, change in the ambient temperature, at a preset time to recalibrate the balance.

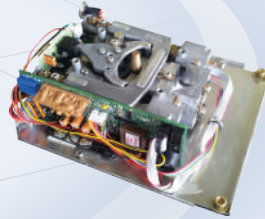

CONTECH[®]
Instruments Ltd.
AN ISO 9001, 14001 & 45001 COMPANY

Where innovation never stops.....

301, Punit Indl. Premises, Turbhe Naka, Navi Mumbai - 400 705.

Tel. : +91 22 6139 3000 (100 lines)

E-mail : sales@contechindia.in

Website: www.contechindia.com

FACTORY / NABL LAB



Plot No. EL-221 TTC Indl. Area, MIDC (Electronic Zone),
Mhape, Navi Mumbai - 400 710.